

### **Giantin Polyclonal Antibody**

Purified Rabbit Polyclonal Antibody (Pab) Catalog # AP55144

## **Specification**

## **Giantin Polyclonal Antibody - Product Information**

Application FC, IHC-P
Primary Accession Q14789
Reactivity Rat, Dog
Host Rabbit
Clonality Polyclonal
Calculated MW 376019

# **Giantin Polyclonal Antibody - Additional Information**

#### **Gene ID 2804**

#### **Other Names**

Golgin subfamily B member 1, 372 kDa Golgi complex-associated protein, GCP372, Giantin, Macrogolgin, GOLGB1

#### **Format**

0.01M TBS(pH7.4), 0.09% (W/V) sodium azide and 50% Glyce

#### Storage

Store at -20  $^{\circ}$ C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4  $^{\circ}$ C.

### **Giantin Polyclonal Antibody - Protein Information**

# Name GOLGB1

#### **Function**

May participate in forming intercisternal cross-bridges of the Golgi complex.

### **Cellular Location**

Golgi apparatus membrane; Single-pass type I membrane protein

## **Giantin Polyclonal Antibody - Protocols**


Provided below are standard protocols that you may find useful for product applications.

- Western Blot
- Blocking Peptides
- Dot Blot
- Immunohistochemistry




- Immunofluorescence
- <u>Immunoprecipitation</u>
- Flow Cytomety
- Cell Culture

# Giantin Polyclonal Antibody - Images

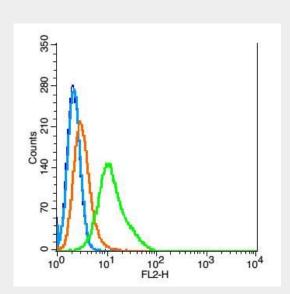

Blank control: Hela(blue), the cells were fixed with 2% paraformaldehyde (10 min) and then permeabilized with ice-cold 90% methanol for 30 min on ice.

Isotype Control Antibody: Rabbit IgG(orange); Secondary Antibody: Goat anti-rabbit IgG-PE(white blue), Dilution: 1:200 in 1 X PBS containing 0.5% BSA; Primary Antibody Dilution: 1 µg in 100 µL1X PBS containing 0.5% BSA(green).

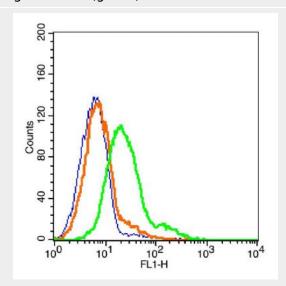


Blank control: Lovo Cells(blue).


Primary Antibody: Rabbit Anti-Giantin/AF488 Conjugated antibody (bs-13356R-AF488), Dilution: 5 µg in 100 µL 1X PBS containing 0.5% BSA;

Isotype Control Antibody: Rabbit IgG/AF488(orange) ,used under the same conditions.




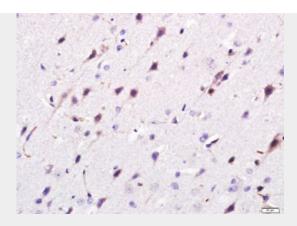

Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Giantin) Polyclonal Antibody, Unconjugated (bs-13356R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.

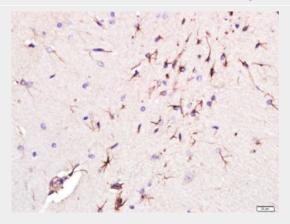


Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Giantin) Polyclonal Antibody, Unconjugated (bs-13356R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.

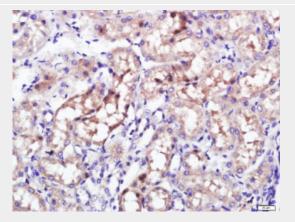


Paraformaldehyde-fixed, paraffin embedded (rat kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Giantin) Polyclonal Antibody, Unconjugated (bs-13356R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.



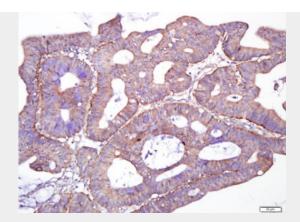

Paraformaldehyde-fixed, paraffin embedded (human cervix cancer); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Giantin) Polyclonal Antibody, Unconjugated (bs-13356R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.

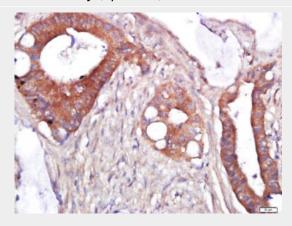

Paraformaldehyde-fixed, paraffin embedded (human cervix cancer); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Giantin) Polyclonal Antibody, Unconjugated (bs-13356R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.