

# **KIAA1755 Polyclonal Antibody**

Purified Rabbit Polyclonal Antibody (Pab)
Catalog # AP56522

# **Specification**

#### **KIAA1755 Polyclonal Antibody - Product Information**

Application IHC-P
Primary Accession O5JYT7
Reactivity Rat
Host Rabbit
Clonality Polyclonal
Calculated MW 130846

## **KIAA1755 Polyclonal Antibody - Additional Information**

**Gene ID 85449** 

#### **Other Names**

Uncharacterized protein KIAA1755, KIAA1755

#### **Format**

0.01M TBS(pH7.4), 0.09% (W/V) sodium azide and 50% Glyce

#### Storage

Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.

## **KIAA1755 Polyclonal Antibody - Protein Information**

Name KIAA1755


# **KIAA1755 Polyclonal Antibody - Protocols**

Provided below are standard protocols that you may find useful for product applications.


- Western Blot
- Blocking Peptides
- Dot Blot
- <u>Immunohistochemistry</u>
- Immunofluorescence
- <u>Immunoprecipitation</u>
- Flow Cytomety
- Cell Culture

# KIAA1755 Polyclonal Antibody - Images



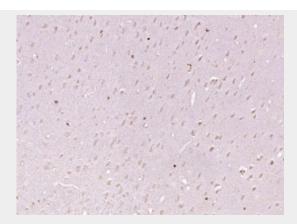

Paraformaldehyde-fixed, paraffin embedded (Rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (KIAA1755) Polyclonal Antibody, Unconjugated (bs-17029R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

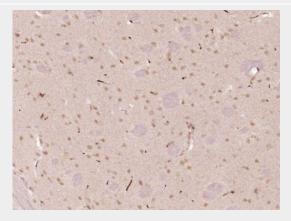

Paraformaldehyde-fixed, paraffin embedded (Mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (KIAA1755) Polyclonal Antibody, Unconjugated (bs-17029R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.