

DRD2 Rabbit pAb

DRD2 Rabbit pAb Catalog # AP94090

Specification

DRD2 Rabbit pAb - Product Information

Application IHC-P
Reactivity Human
Host Rabbit
Clonality Polyclonal

DRD2 Rabbit pAb - Additional Information

Format

0.01M TBS(pH7.4), 0.09% (W/V) sodium azide and 50% Glyce

Storage

Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.

DRD2 Rabbit pAb - Protein Information


DRD2 Rabbit pAb - Protocols

Provided below are standard protocols that you may find useful for product applications.


- Western Blot
- Blocking Peptides
- Dot Blot
- Immunohistochemistry
- Immunofluorescence
- <u>Immunoprecipitation</u>
- Flow Cytomety
- Cell Culture

DRD2 Rabbit pAb - Images




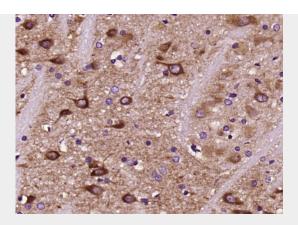

Paraformaldehyde-fixed, paraffin embedded (Rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Dopamine D2 Receptor) Polyclonal Antibody, Unconjugated (AP94090) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

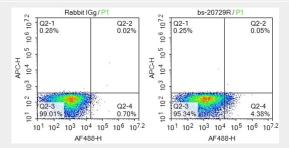


Blank control:Mouse brain. Primary Antibody (green line): Rabbit Anti-Dopamine Receptor D2 antibody (AP94090) Dilution: 2 μ g /10^6 cells; Isotype Control Antibody (orange line): Rabbit IgG . Secondary Antibody: Goat anti-rabbit IgG-AF488 Dilution: 1 μ g /test. Protocol The cells were incubated in 5%BSA to block non-specific protein-protein interactions for 30 min at room temperature. Cells stained with Primary Antibody for 30 min at room temperature. The secondary antibody used for 40 min at room temperature. Acquisition of 20,000 events was performed.

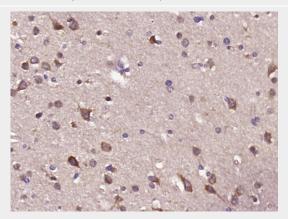

Paraformaldehyde-fixed, paraffin embedded (Human glioma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Dopamine D2 Receptor) Polyclonal Antibody, Unconjugated (AP94090) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

DRD2 Rabbit pAb - Background

This product as supplied is intended for research use only, not for use in human, therapeutic or

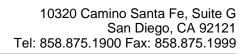



diagnostic applications.