

#### liver FABP Rabbit pAb

liver FABP Rabbit pAb Catalog # AP94444

## **Specification**

## liver FABP Rabbit pAb - Product Information

Application IHC-P, IHC-F, IF
Reactivity Chicken
Host Rabbit
Clonality Polyclonal

## liver FABP Rabbit pAb - Additional Information

#### **Dilution**

<span class ="dilution\_IHC-P">IHC-P~~N/A</span><br \> <span class = "dilution\_IHC-F">IHC-F~~N/A</span><br \> <span class = "dilution\_IF">IF~~1:50~200</span>

#### **Format**

0.01M TBS(pH7.4), 0.09% (W/V) sodium azide and 50% Glyce

## **Storage**

Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.

## liver FABP Rabbit pAb - Protein Information


## liver FABP Rabbit pAb - Protocols

Provided below are standard protocols that you may find useful for product applications.


- Western Blot
- Blocking Peptides
- Dot Blot
- <u>Immunohistochemistry</u>
- Immunofluorescence
- Immunoprecipitation
- Flow Cytomety
- Cell Culture

# liver FABP Rabbit pAb - Images




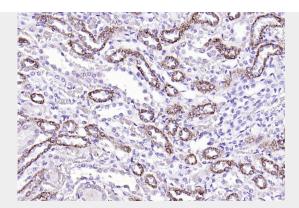

Paraformaldehyde-fixed, paraffin embedded (rat kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (liver FABP) Polyclonal Antibody, Unconjugated (AP94444) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining

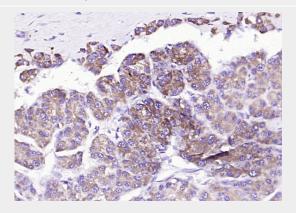


Paraformaldehyde-fixed, paraffin embedded (human liver carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (liver FABP) Polyclonal Antibody, Unconjugated (AP94444) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

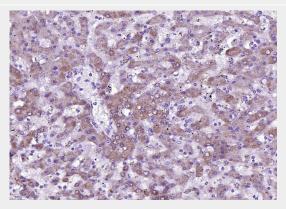


Paraformaldehyde-fixed, paraffin embedded (human liver); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (liver FABP) Polyclonal Antibody, Unconjugated (AP94444) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.




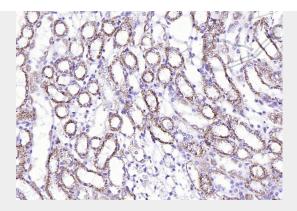

Paraformaldehyde-fixed, paraffin embedded (mouse kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (liver FABP) Polyclonal Antibody, Unconjugated (AP94444) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining

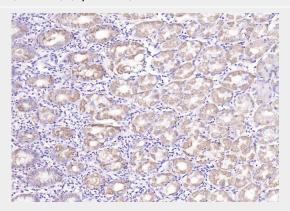

Paraformaldehyde-fixed, paraffin embedded (Human stomach); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (liver FABP) Polyclonal Antibody, Unconjugated (AP94444) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

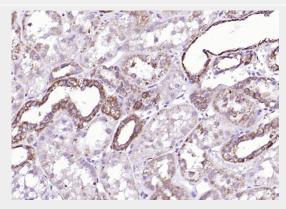

Paraformaldehyde-fixed, paraffin embedded (Human kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (liver FABP) Polyclonal Antibody, Unconjugated (AP94444) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

liver FABP Rabbit pAb - Background





This product as supplied is intended for research use only, not for use in human, therapeutic or diagnostic applications.