

Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody

Rabbit Polyclonal, Peroxidase (Horseradish)
Catalog # ASR1946

Specification

Physical State

Immunogen

Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Product Information

Description Anti-RAT IgG (H&L) (RABBIT) Antibody

Peroxidase Conjugated

Host Rabbit

Conjugate Peroxidase (Horseradish)

Target Species Rat

Clonality Polyclonal Application ,1,10,15,

Application Note ELISA 1:170,000; Western Blot

1:1,000-1:5,000;Immunochemistry

1:500-1:2,500 Lyophilized

Host Isotype IgG

Target Isotype IgG (H&L)

Buffer 0.02 M Potassium Phosphate, 0.15 M

Sodium Chloride, pH 7.2 Rat IgG whole molecule

Reconstitution Volume 1.0 mL

Reconstitution Buffer Restore with deionized water (or

equivalent)

Stabilizer 10 mg/mL Bovine Serum Albumin (BSA) -

Immunoglobulin and Protease free

Preservative 0.01% (w/v) Gentamicin Sulfate. Do NOT

add Sodium Azide!

Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Additional Information

Shipping Condition

Ambient

Purity

Secondary Antibody Conjugated was prepared from monospecific antiserum by immunoaffinity chromatography using Rat IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Peroxidase, anti-Rabbit Serum, Rat IgG and Rat Serum.

Storage Condition

Store secondary antibody at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.

Precautions Note

This product is for research use only and is not intended for therapeutic or diagnostic applications.




Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Protein Information

Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Protocols

Provided below are standard protocols that you may find useful for product applications.

- Western Blot
- Blocking Peptides
- Dot Blot
- <u>Immunohistochemistry</u>
- Immunofluorescence
- Immunoprecipitation
- Flow Cytomety
- Cell Culture

Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Images

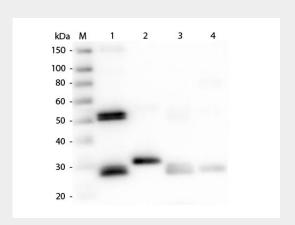

Western Blot of Anti-Rat IgG (H&L) (RABBIT) Antibody . Lane M: 3 μ l Molecular Ladder. Lane 1: Rat IgG whole molecule . Lane 2: Rat IgG F(c) Fragment . Lane 3: Rat IgG F(ab) Fragment . Lane 4: Rat IgM Whole Molecule . All samples were reduced. Load: 50 ng per lane. Block: MB-070 for 30 min at RT. Primary Antibody: Anti-Rat IgG (H&L) (RABBIT) Antibody 1:1,000 for 60 min at RT. Secondary Antibody: Anti-Rabbit IgG (GOAT) Peroxidase Conjugated Antibody 1:40,000 in MB-070 for 30 min at RT. Predicted/Obsevered Size: 25 and 55 kDa for Rat IgG, 25 kDa for F(c) and F(ab), 78 and 25 kDa for IgM. Rat F(c) migrates slightly higher.

Anti-RAT IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Background

HRP secondary antibody conjugated Rat IgG is Suitable for western blot, ELISA, immunoperoxidase electron microscopy and immunohistochemistry as well as other peroxidase-antibody based enzymatic assays requiring lot-to-lot consistency.