

Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody

Goat Polyclonal, Peroxidase (Horseradish)
Catalog # ASR2962

Specification

Host

Physical State

Host Isotype

Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody - Product Information

Description Anti-RAT IgG (H&L) (GOAT) Antibody

Peroxidase Conjugated (Min X Bv Ch Gt GP Ham Hs Hu Ms Rb & Sh Serum Proteins)

Goat

Conjugate Peroxidase (Horseradish)

Target Species
Reactivity
Clonality
Application
Reit Rat
Peroxidase
Rat
Rat
Polyclonal
WB, E, IC

Application Note ELISA 1:100,000; Western Blot

1:2,000-1:10,000; Immunochemistry

1:500-1:1,500 Lyophilized

IqG

Target Isotype IgG (H&L)

Buffer 0.02 M Potassium Phosphate, 0.15 M

Sodium Chloride, pH 7.2

Immunogen Anti-Rat IgG whole molecule was produced

by repeated immunization with Rat IgG

whole molecule in goat.

Reconstitution Volume 1.0 mL

Reconstitution Buffer Restore with deionized water (or

equivalent)

Stabilizer 10 mg/mL Bovine Serum Albumin (BSA) -

Immunoglobulin and Protease free

Preservative 0.01% (w/v) Gentamicin Sulfate. Do NOT

add Sodium Azide!

Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody - Additional Information

Shipping Condition

Ambient

Purity

This product was prepared from monospecific antiserum by immunoaffinity chromatography using Rat IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Peroxidase, anti-Goat Serum, Rat IgG and Rat Serum. No reaction was observed against Bovine, Chicken, Goat, Guinea Pig, Hamster, Horse, Human, Mouse, Rabbit and Sheep Serum Proteins.

Storage Condition



Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.

Precautions Note

This product is for research use only and is not intended for therapeutic or diagnostic applications.


Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody - Protein Information

Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody - Protocols

Provided below are standard protocols that you may find useful for product applications.

- Western Blot
- Blocking Peptides
- Dot Blot
- Immunohistochemistry
- Immunofluorescence
- Immunoprecipitation
- Flow Cytomety
- Cell Culture

Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody - Images

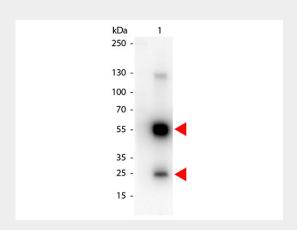

Western Blot of Peroxidase conjugated Goat anti-Rat IgG antibody. Lane 1: Rat IgG. Lane 2: none. Load: 50 ng per lane. Primary antibody: none. Secondary antibody: Peroxidase rat secondary antibody at 1:1,000 for 60 min at RT. Block: MB-070 for 30 min at RT. Predicted/Observed size: 55 kDa, 28 kDa for Rat IgG. Other band(s): none.





Western Blot of Anti-Rat IgG (H&L) (GOAT) Antibody (Min X Bv Ch Gt GP Ham Hs Hu Ms Rb & Sh Serum Proteins) . Lane M: 3 μ l Molecular Ladder. Lane 1: Rat IgG whole molecule . Lane 2: Rat IgG F(c) Fragment . Lane 3: Rat IgG F(ab) Fragment . Lane 4: Rat IgM Whole Molecule . Lane 5: Rat Serum . All samples were reduced. Load: 50 ng per Iane. Block: MB-070 for 30 min at RT. Primary Antibody: Anti-Rat IgG (H&L) (GOAT) Antibody (Min X Bv Ch Gt GP Ham Hs Hu Ms Rb & Sh Serum Proteins) 1:1,000 for 60 min at RT. Secondary Antibody: Anti-Goat IgG (DONKEY) Peroxidase Conjugated Antibody 1:40,000 in MB-070 for 30 min at RT. Predicted/Obsevered Size: 25 and 55 kDa for Rat IgG and Serum, 25 kDa for F(c) and F(ab), 78 and 25 kDa for IgM. Rat F(c) migrates slightly higher.

Anti-Rat IgG (H&L) (Peroxidase Conjugated) Pre-Adsorbed Secondary Antibody - Background

Anti-Rat IgG whole molecule antibody generated in goat detects specifically Rat IgG whole molecule. This secondary antibody anti-Rat is ideal for investigators who routinely perform ELISA, Sandwich ELISA, titration assays, western-blot, immunoprecipitation and more generally immunoassays.